




**Interface / Receiver**  
(Ordering Code 6070100)



All dimensions in mm

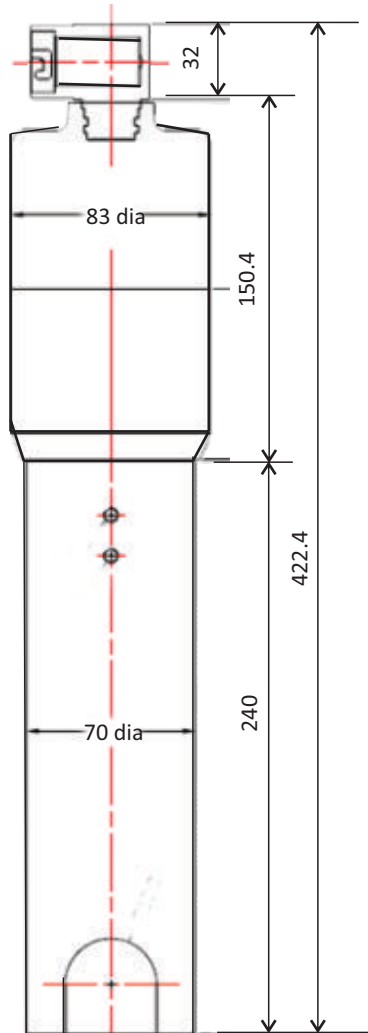

Kink-resistance transfer hose with a pair of 30 mm conical (ISO) Male & Female connector.  
Male connector at Interface end & Female connector at anaesthesia workstation exhaust port.  
(Ordering Code 6070105)



3m BS Compatible AGSS Evacuation hose assembly (5m on request)  
(Ordering Code 6071100)

## Description

Scavenging is the collection and removal of vented anaesthetic gases from the Operating Room and this receiver fulfills the vital function of protecting the patient from excessive positive or negative pressures. Positive-pressure relief is mandatory to vent excess gas in case of occlusion distal to a closed interface. It is the most important component in an active system for protection from excessive positive system pressure and reduction of waste anaesthetic gases (WAGs) in an operating room environment to avoid occupational exposure.

The Scavenging Interface / receiver is placed between breathing and disposal system i.e. Transfer hose which transports WAGs from the exhaust valve of an anaesthetic ventilator or anaesthetic breathing system and Evacuation hose (Disposal hose) connected to the AGSS Outlet for eventual disposal into the open atmosphere at a safe distance.



**Interfacing / Receiving Unit Set-up**

## Ordering and Packing Information

**Ordering Code:** 6070101

Interfacing / Receiving Unit comprising of:

- a. Receiver
- b. Kink-resistance transfer hose
- c. A pair of 30 mm conical (ISO) Male & Female connector
- d. 3m Evacuation / disposal hose with crimped probe & Connector (5m on request)



**AKTIV TECHNOLOGIES**

**Corporate Office :** Plot No. 691, Sector 69, IMT, Faridabad, Haryana 121004, INDIA

Tel: +91-7291838253

Email: [aktivtechnologies@gmail.com](mailto:aktivtechnologies@gmail.com), [sanjay.aktiv@gmail.com](mailto:sanjay.aktiv@gmail.com)

Website : [www.aktivtechnologies.com](http://www.aktivtechnologies.com)



Consult Instructions  
for Use

LT / 607A Rev 00 12/21

**Disclaimer:** Due to our policy of continuous product improvement, Aktiv Technologies reserve the right to alter dimensions and/or specification without notice although every effort will be made to advise of any such modifications.